



Mechanical Baby Scales

- Time-tested, reliable weigh beam accuracy
- No electricity required – ideal for remote clinics
- Weigh beam displays readings on both sides
- Baked-on, powder painted steel base and column
- Comfortable, polystyrene weighing tray
- Built-in measuring tape (up to 21 in / 53 cm)

DETECTO
450 | 451

DETECTO's model 450 baby scale shown here with 130 lb x 1 oz capacity and built-in measuring tape.



DISPLAYS FROM EITHER SIDE

The weigh beam reading displays from either side, front and back, providing superior ease of use for nurses.



TIME-TESTED WEIGHBEAM

DETECTO's die-cast mechanical weigh beam is world-renowned for its high accuracy and long-term dependability over decades of use.



BUILT-IN MEASURING TAPE

The easy-to-read measuring tape allows users to measure both weight and height simultaneously (0 - 21 in/53 cm).

Mechanical Baby Scales

2 models to choose from • Durable, clinical quality

DETECTO's 450 series mechanical baby scales provide the ultimate in safety for pediatric patients while also utilizing the time-tested, reliable DETECTO weigh beam for highly-accurate measurement. They are available in capacities up to 130 pounds and 65 kilograms with your choice of pounds or kilograms. There is no electricity required, making them ideal for remote clinics and pediatrician offices.

The easy-to-clean polystyrene tray cradles babies securely on these dependable pediatric scales. The heavy-duty steel base provides patient assurance while weighing. A built-in measuring tape allows the healthcare professional to measure both weight and height simultaneously.


The die-cast weigh beam is positioned high above the scale base to provide for unobstructed viewing of the baby's weight measurement. The baked-on, white powder paint offers a long-lasting dependable finish. For fast mechanical measurement, DETECTO is the most-trusted name brand in clinical weighing.



SPECIFICATIONS

Capacities	450: 130 lb x 1 oz 451: 65 kg x 20 g
Overall Dimensions	22 in W x 21 in D x 27 in H / 56 cm W x 53 cm D x 69 cm H
Tray Dimensions	22 in W x 14.75 in D x 4.5 in H / 56 cm W x 37 cm D x 11 cm H
Construction	Base and Column: White baked-on, powder-painted steel Weighing Tray: White polystyrene
Measuring Tape	0 - 21 in / 0 - 53 cm
Shipping Weight	33 lb / 15 kg

Detecto Scale reserves the right to improve, enhance, or modify features and specifications without prior notice.

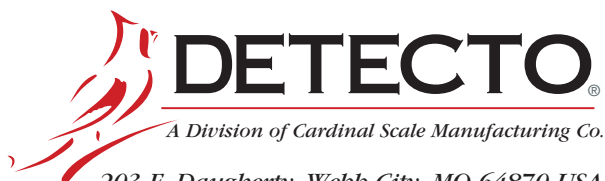

DETECTO[®]

A Division of Cardinal Scale Manufacturing Co.

203 E. Daugherty, Webb City, MO 64870 USA

Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-5001

www.Detecto.com

Mobile: m.detecto.com

SOLD BY: