

model **7550**
Wall-Mount Fold-up
WHEELCHAIR
SCALES



WiFi
Optional

SPACIOUS PLATFORM

32.3-in W x 33-in D /
82-cm W x 84-cm D
Weighing Platform
Provides Ample Room

Detecto
Model 7550
Wheelchair Scale
with Easy-Access Ramps

- 1,000 lb / 300 kg Capacity
- Optional Wi-Fi for EMR/EHR
- Sturdy Steel Construction




1,000-lb/300-kg CAPACITY

Ideal for Bariatric Weighing



Detecto's 7550 Wheelchair Scale folds up compact against the wall to conserve floor space when not in use.

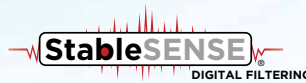

- Compact Vertical Fold-Up Platform Conserves Space
- Pushbutton and Keypad Tare for Wheelchair Weights
- Accommodates Some Sizes of Power Wheelchairs
- May be Used on Carpet or Hard Floor Surfaces
- RS232 Serial and USB Ports Standard
- Body Mass Index Calculation
- P600 Tape Printer Included with Model 7550-PR-AC



Detecto's model 7550 wheelchair scale features a unique fold-up storage design that saves space when not in use.

The 7550 easily accommodates most wheelchairs, and yet requires very little floor space when folded up. The slim 6.6-inch / 17-cm profile takes up minimal space when installed in a medical facility hallway, clinical office, or patient room. When you're ready to weigh, simply lower the platform and the gas spring lifting mechanism will gently lower it into place. The patient may be easily rolled onto the scale, and a platform lip will hold the wheelchair securely in place for weighing.




Slim 6.6-inch / 17-cm Compact Profile

MOUNTS SECURELY TO A WALL FOR STORAGE

Compact Storage

After use, the 7550 may be folded up vertically to conserve space. Patients may be weighed in a wheelchair, chair, or standing.

Optional Wi-Fi for EMR/EHR 

Tilt-and-Swivel Rotating Indicator



Sturdy Steel Construction for Durability

Gas Spring Platform Lifting Mechanism

Textured, Non-Skid Mat on Ramp and Platform

Two-Way Easy Access Ramps

Platform Lip to Hold Wheelchairs in Place on Scale

EFFICIENT MEASUREMENT



The 7550's MedVue® weight indicator is the clinical industry's most advanced, featuring EMR/EHR-ready serial, USB, and optional Wi-Fi connectivity, BMI calculation, yellow "quick keys" for basic weighing, and a multi-line LCD display. The scale is powered by 6 C batteries (not included) or optional AC adapter.



model **7550** Wall-Mount Fold-up **WHEELCHAIR SCALES**



DETECTO is the largest medical scale manufacturer in the United States providing clinical-grade accuracy and dependability for over 113 years.



The 7550's 1.5-inch/38-mm high platform is very low to the ground for easy patient access. This makes the process simple and minimal effort required in order to roll a patient on and off the weighing platform.

Model	7550	7550-AC	7550-PR-AC
Description	AC adapter and printer not included	AC adapter included	AC adapter and P600 tape printer with 8545-B098-0A connecting cable included
Capacity	1,000 lb x 0.2 lb / 300 kg x 0.1 kg		
Platform Size (Flat Portion)	32.3 in W x 33 in D x 1.5 in H / 81.9 cm W x 83.9 cm D x 3.8 cm H		
Overall Dimensions	Width: 50.3 in / 127.6 cm Depth: 39.3 in / 99.7 cm (with platform extended) 6.6 in / 16.8 cm (without platform extended) Height: 48.8 in / 123.8 cm (with indicator) 41.1 in / 104.4 cm (without indicator)		
Weight Indicator Dimensions	9.8 inch W x 6.8 inch H x 1.8 inch D / 248.9 mm W x 172.7 mm H x 45.7 mm D		
LCD Display Size	3.5 inch W x 3 inch H / 89 mm W x 76 mm H backlit blue LCD		
Display	Weight: 5 digit, seven segment, 0.875 inch / 22.2 mm high LCD Height: 4 digit, seven segment, 0.5 inch / 12.7 mm high LCD Body Mass Index: 3 digit, seven segment, 0.5 inch / 12.7 mm high LCD ID: 14 digit, 14 segment, 0.35 inch / 8.9 mm high LCD		
Resolution	5,000 divisions, maximum		
Keypad	Sealed-membrane type with 19 color-coded keys		
Connectivity Ports	RS232 serial and USB (standard), wired and Wi-Fi Ethernet (optional)		
Power	6 C cell Alkaline, Ni-Cad or NiMH batteries (not included) OR 12 VDC 1.25A wall plug-in AC power adapter, Detecto part number MV1PWR (included with -AC models listed above)		
Battery Life	7,560 weighments (with backlight fully ON) / 10,620 weighments (with backlight fully OFF)		
Weighing Mechanism	Electronic load cells		
Lifting Mechanism	Gas springs		
Operating Temperature	14° to 104° F (-10° to +40° C)		
Construction	Steel		
UPC Code	809161190908	809161316803	809161316902
Net Weight	180 lb / 82 kg	181 lb / 82 kg	184 lb / 83 kg
Shipping Weight	238 lb / 108 kg	239 lb / 109 kg	242 lb / 110 kg

Detecto Scale reserves the right to improve, enhance, or modify features and specifications without prior notice.

 **DETECTO**
A Division of Cardinal Scale Manufacturing Co.
102 E. Daugberty, Webb City, MO 64870 USA
Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-2153
www.Detecto.com

Mobile: m.detecto.com

SOLD BY:

