

PD series

PROFESSIONAL DOCTOR SCALES


MEDICAL GRADE ACCURACY: WEIGHT • HEIGHT • BMI



Detecto Model PD300M



Detecto Model PD100M
Low-Profile Physician Scale
with Body Mass Index

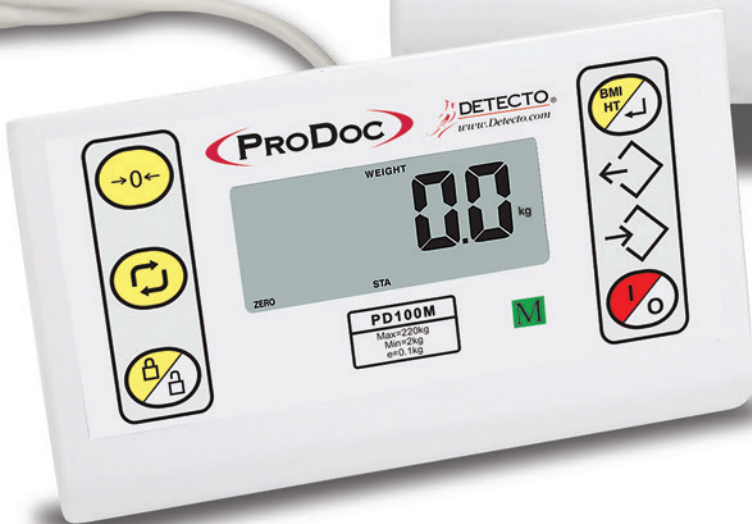

Detecto Model PD100M Low-Profile Physician Scale with Remote Display



Model PD100M

Low-Profile Floor Scale with Remote Display

The PD100M's remote display may be positioned on desks, counters, or even mounted on the wall (shown below) up to 220 cm from the scale. The LCD provides easy-to-view 25 mm readouts. Include the optional DHR wall-mount digital

height rod and model GBSS18-WM stainless steel handrails to create an all-digital clinical measurement station.



Optional model PRODOC-CASE carrying case for PD100M

Model PD200M

Space-Saving Physician Scale

The PD200M features a low-profile, angled display for easy measurement viewing. The PD100M and PD200M provide BMI (Body Mass Index) calculation by conveniently selecting your height with easy-to-use scroll keys while weighing.



Model PD200M Space-Saving Physician Scale



Model DHR Digital Height Rod with Wall-Mount Brackets

Model GBSS18-WM Stainless Steel Wall-Mount Handrails

Model PD100M Low-Profile Physician Scale with Body Mass Index



Optional P50 printer may be mounted to any of the PD300M series columns

09:28, 07/30/12
 BMI: 21.5
 Height: 175.5 cm
 Weight: 68.5 kg

EYE-LEVEL MEASUREMENT CONVENIENCE

The optional P50 direct thermal printer may be mounted on any of the ProDoc® PD300M series eye-level physician scales to print Time, Date, Height, Weight, and Body Mass Index. The P50 features two-button operation and RS-232 serial communication with the scale (cable included with printer).




Model PD300M
 Eye-Level Physician Scale
 with Body Mass Index



Digital
 Height Rod
 Measures from
 110 - 200 cm

Model PD300MDHR
 Eye-Level Physician Scale
 with BMI and Digital Height Rod

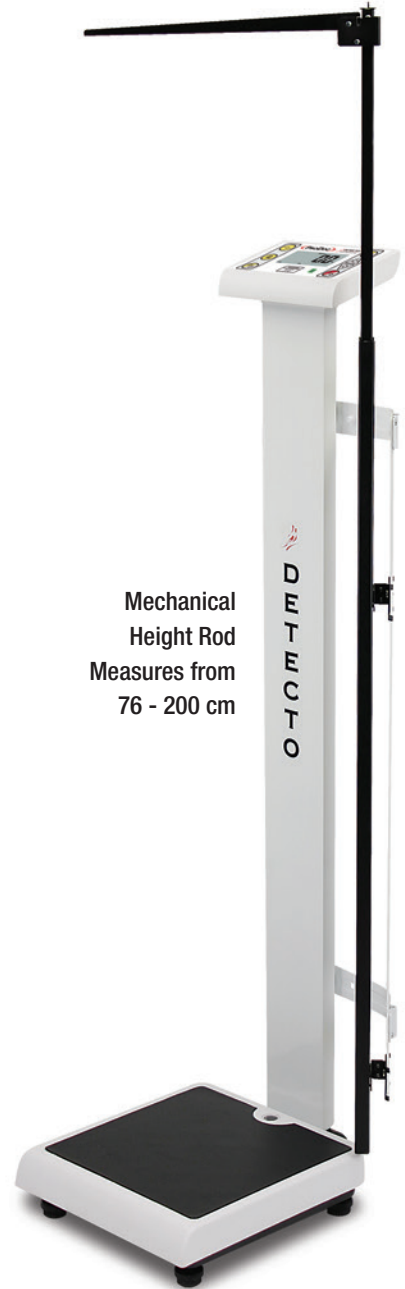

Mechanical
 Height Rod
 Measures from
 76 - 200 cm

Model PD300MMHR
 Eye-Level Physician Scale
 with BMI and Mechanical Height Rod



PROFESSIONAL DOCTOR SCALES

Level Indicator

Adjustable Locking Feet



Built-in wheels for fast and easy transport



Deteco's ProDoc® Series Digital Physician Scales with Body Mass Index


offer medical-grade accuracy and durability from the brand you trust in clinical weighing at an economical price. The PD300MDHR features time-saving digital weight and height measurement, so the patient simply steps on the scale, raises the height rod into position, and presses a key to automatically calculate their BMI. The PD300M series features comfort-height-level display of measurements for convenient viewing.

Multiple configurations to fit your weighing application

- Capacity: 250 kg - PD300M series; 220 kg - PD100M/200M
- 4-user memories and recall function
- 6 AA batteries (not included) or optional AC adapter
- Solid, stable base for patient assurance and accuracy
- Height measurement (available on models equipped with height rod)
- Body Mass Index (BMI) calculation
- Metric conversion (Weight and Height)
- Selectable Sleep Mode (up to 9 minutes)
- Selectable Auto Shut-off Feature (up to 9 minutes)
- Serial and USB ports on PD300M series
- Available with or without height rod
- Large 25-mm high LCD readouts




PD300M Series Keypad Features

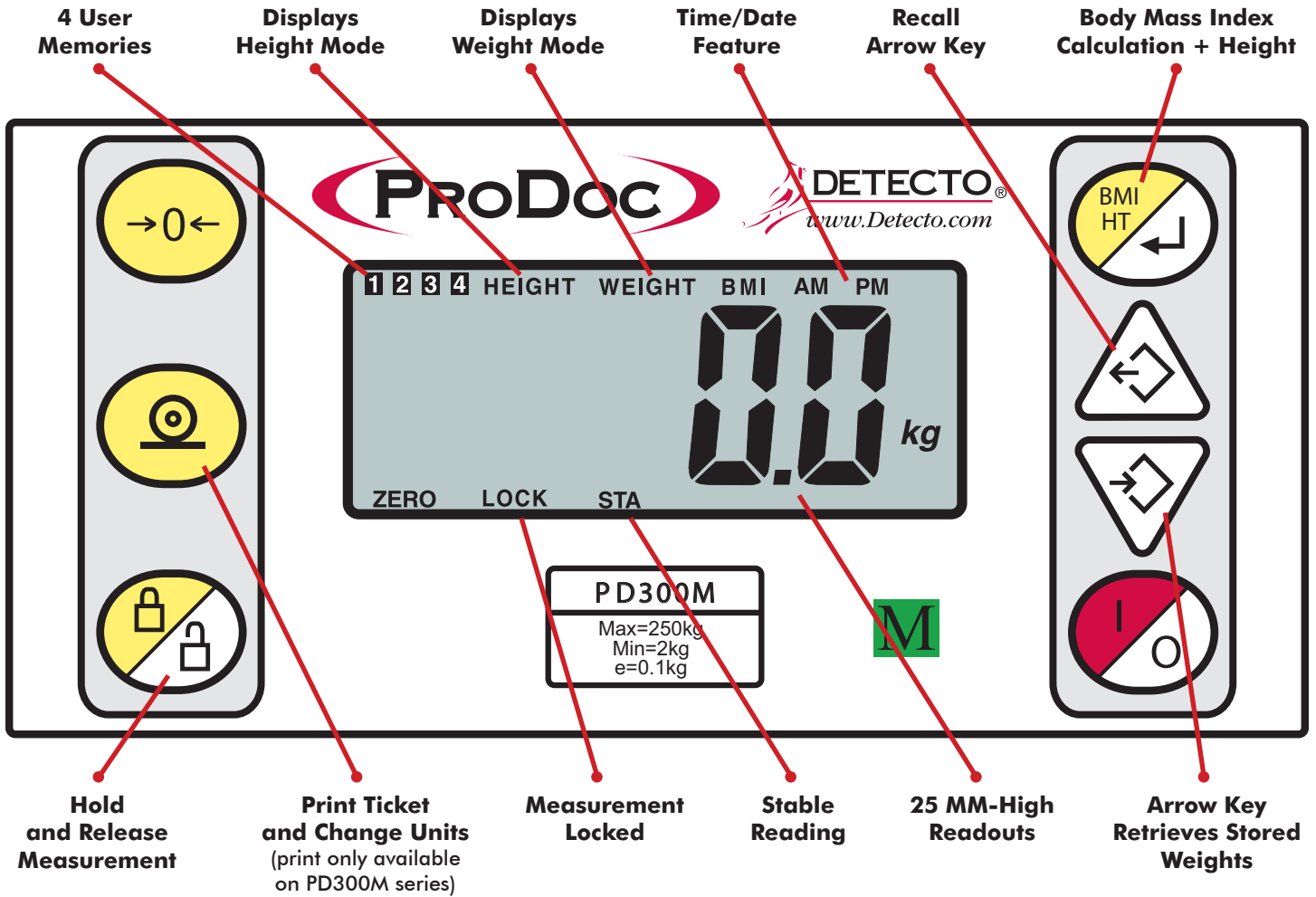

Features USB and Serial Connectivity for Efficient Patient Records

Back View of Indicator (PD300M Series Only Features USB and Serial Ports / Not Available on PD100M-200M)

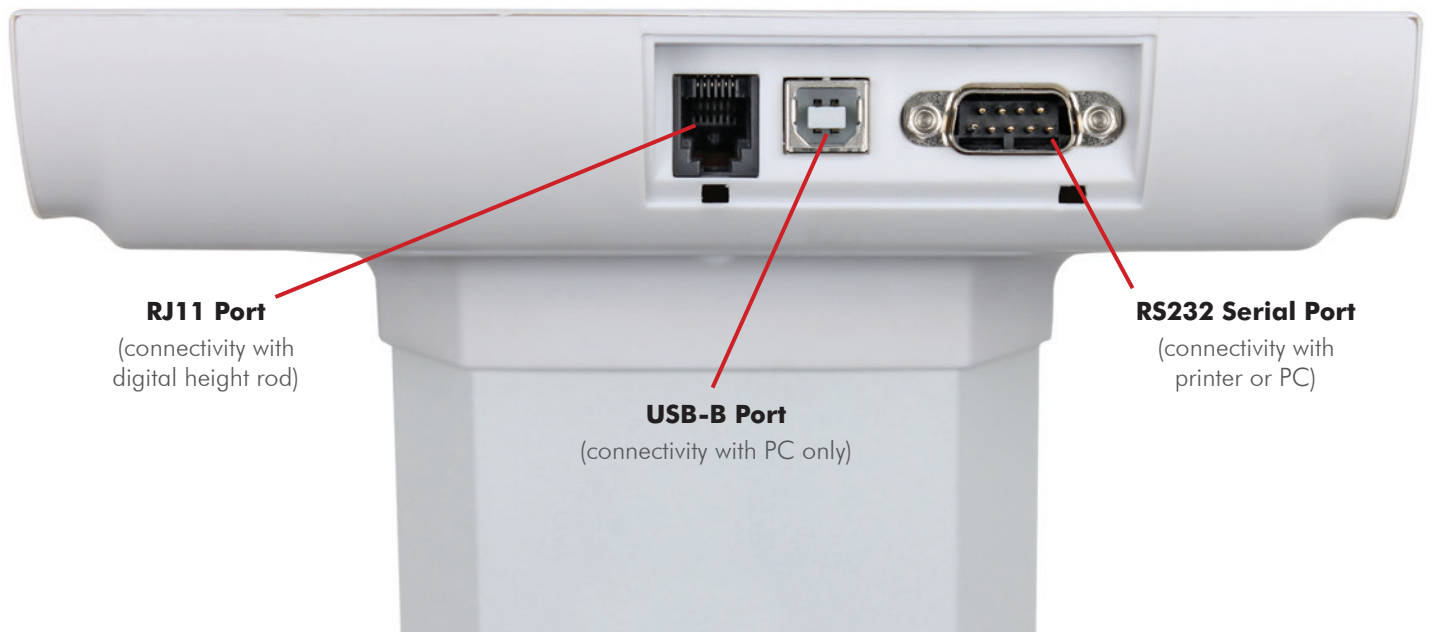



PD series
**PROFESSIONAL
DOCTOR SCALES**

PD SERIES SCALE SPECIFICATIONS

Model Number	PD100M (Floor Scale with Remote Display) PD200M (Floor Scale) PD300M (Eye-Level Scale) PD300MDHR (Eye-Level Scale with Digital Height Rod) PD300MMHR (Eye-Level Scale with Mechanical Height Rod)
Capacity PD100M/PD200M PD300M	220 kg x 0.1 kg 0-100 kg x 50 g 100-250 kg x 100 g
Weight Units	Kilograms
Power Requirements	Six (6) "AA" size Alkaline batteries (not included) OR an optional Medical device 9V AC/DC wall plug-in adapter. Cardinal part number: UK plug (PD-UKPLUG) or EU plug (PD-EUPLUG).
Display	Five digit, seven segment, 25 mm high LCD
Operating Temperature	+10 to +40° C
Function Keys	ON/OFF, BMI HEIGHT/ENTER, LOCK/RELEASE, UNITS (PRINT/UNITS on PD300M), ZERO, RECALL and STORE
Dimensions PD100M Display PD100M Base PD200M PD300M PD300MDHR PD300MMHR	134 mm x 200 mm x 54 mm 85 mm x 318 mm x 330 mm 222 mm x 320 mm x 430 mm 1,219 mm x 320 mm x 433 mm 1,219 mm x 354 mm x 433 mm 1,219 mm x 354 mm x 433 mm

P50 PRINTER SPECIFICATIONS (for use with PD300M series)

Printing Type	Direct Thermal Dot Line
Number of Dots/Line	384 dots
Resolution	8 dots/mm W x 16 dots/mm H
Max. Printing Speed	40 mm/second
Max. Paper Roll Diameter	50 mm
Printing Width	48 mm
Operating Temperature	20 to 50° C
Power Supply	AC Adapter 100 - 240VAC input, 12V [2A] output
Connectivity	RS-232 communication with scale (cable included)
Dimensions	110 mm x 110 mm x 107 mm

Detecto Scale reserves the right to improve, enhance, or modify features and specifications without prior notice.



DETECTO[®]
A Division of Cardinal Scale Manufacturing Co.
203 E. Daugherty, Webb City, MO 64870 USA
Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-5001
www.Detecto.com

Sold By: