


# STRETCHER SCALE®



*Detecto  
Model  
8550*

## MODELS 8500/8550

- Spacious platform accommodates all common stretcher types
- Two integral wheels and guide handles for mobility
- High-capacity digital scale for stretchers and wheelchairs
- Available with or without column for weight indicator

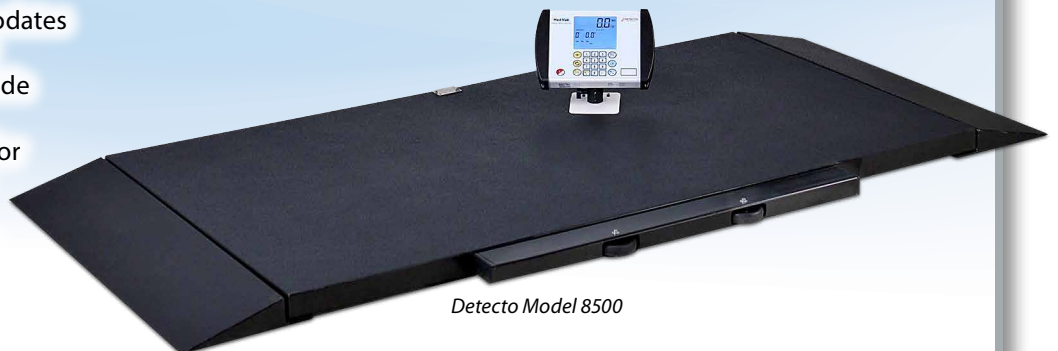

*Detecto Model 8500*

# PRECISION WEIGHING

## FOR PATIENTS IN STRETCHERS GURNEYS OR WHEELCHAIRS



Detecto  
Model  
8500



Detecto  
Model 8550

- 800 lb x 0.2 lb / 360 kg x 0.1 kg capacity
- Spacious 60 in W x 32.5 in D / 1,524 mm W x 826 mm D platform size
- Hold/Release key to hold weight after leaving the scale for better patient care
- Sturdy steel construction for optimal durability
- Mobile for use in multiple clinical departments
- Two integrated ramps included for easy access



**DETECTO**  
[www.Detecto.com](http://www.Detecto.com)

DETECTO's Stretcher Scale® is ideal for emergency arrivals into a hospital when time is of the essence for gathering patient weights. The DETECTO models 8500 and 8550 can be used in the ER, dialysis, long-term care, and nursing homes.

The single-pass weighing feature saves time by allowing you to enter the stretcher or wheelchair tare weight using the numeric keypad on the indicator. You may also use the convenient pushbutton tare feature to remove the stretcher or wheelchair weight if you don't know the weight already.

**StableSENSE**  
DIGITAL FILTERING

**EMR  
EHR  
ready**

**WelchAlllyn**  
OPTIONAL

**Wi Fi**  
OPTIONAL

### PORTABLE WITH INTEGRAL WHEELS AND GUIDE HANDLES

DETECTO's Stretcher Scale® is engineered and manufactured at a light enough weight that it is easily portable to lift up and roll by nursing professionals. The model 8500 is 115 lb/ 52 kg net weight and the model 8550 is 130 lb/59 kg. The slim profile scale design fits easily through doorways of nearly any size.





# DETECTO'S STRETCHER SCALE® IS LOADED WITH HELPFUL FEATURES AND TIME-SAVING BENEFITS



## STRETCHER SCALE®



### LIFT AND ROLL

The scale's relatively low net weight, wheels, and guide handles all combine to make it easy to maneuver.



### WEIGH WHEELCHAIRS, TOO

The Stretcher Scale® may be used to weigh wheelchairs as well as stretchers and gurneys.



### LOW PROFILE PLATFORM

DETECTO's Stretcher Scale® has a 2-inch-high, low-profile platform with textured mat surface.



### BRACKET INCLUDED (Model 8500)

The model 8500 comes with a wall or desk mount swivel bracket for the MedVue® weight indicator.



### EASILY-COLLAPSIBLE COLUMN

It's easy to quickly slide the release lever and fold down the 8550's column for transport.




### TWO-WAY QUICK ACCESS RAMPS

The two-way-access ramps rest on the floor making it easy to roll patients onto the scale platform.



### OPTIONAL GUIDE RAILS

An optional yellow guide rail kit is available for the 8500/8550 (8500-GUIDERAILS).

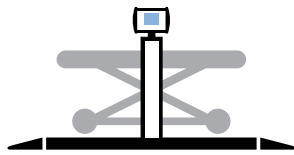

# STRETCHER SCALE®

## DETECTO MODELS 8500/8550



DETECTO model 8550 shown with optional yellow guide rails

### Options

<b>8500-GUIDERAILS</b>	Guide rail kit painted yellow for models 8500/8550 (includes wheel side and handle side)
<b>MV1PWR</b>	Medical-grade AC adapter power supply with US plug
<b>MVWIFI</b>	Wireless option card
<b>MVIP</b>	Wired Ethernet interface option



<b>Capacity</b>	800 lb x 0.2 lb / 360 kg x 0.1 kg	
<b>LCD Display Dimensions</b>	3.5 in W x 3 in H / 89 mm W x 76 mm H backlit blue LCD	
<b>Display</b>	5 digit, 7 segment, 7/8 inch (22 mm) high LCD weight digits	
<b>Resolution</b>	5,000 divisions, max.	
<b>Keypad</b>	19 key, color-coded, sealed membrane	
<b>Functions</b>	Tare, Net/Gross, Hold/Release, Zero, metric conversion, ID storage, Height entry for Body Mass Index calculation, and printer output	
<b>Power</b>	Optional model MV1PWR AC adapter or 6 C batteries (not included)	
<b>Battery Life</b>	7,560 weighments (with backlight fully ON) 10,620 weighments (with backlight fully OFF)	
<b>Connectivity</b>	RS232 serial, USB and optional wired and Wi-Fi Ethernet	
<b>Cable Length (Model 8500 only)</b>	6 ft / 1.8 m coiled cable, 14 ft / 4.2 m fully extended (from scale to indicator)	
<b>Mounting Bracket</b>	350-degree indicator swivel mounting bracket included	
<b>Operating Temperature</b>	14 to 104 °F / -10 to +40 °C	
<b>Construction</b>	Mild steel	
<b>Included</b>	Two ramps (assembly required)	
<b>Overall Dimensions</b>	76 in W x 38 in D x 2 in H / 1,930 mm W x 965 mm D x 51 mm H	76 in W x 38 in D x 44.8 in H / 1,930 mm W x 965 mm D x 1,138 mm H
<b>Remote Indicator Dimensions</b>	9.8 in W x 8.3 in D x 7.8 in H / 249 mm W x 211 in D x 198 mm H (including swivel mounting bracket)	Included in overall dimensions
<b>Platform Dimensions</b>	60 in W x 32.5 in D x 2 in H / 1,524 mm W x 826 mm D x 51 mm H	
<b>Ramp Dimensions</b>	7.5 in W x 32.5 in D x 2 in H / 191 mm W x 826 mm D x 51 mm H	
<b>Optional Guide Rail Dimensions</b>	24 in W x .1 in D x 3.3 in H / 610 mm W x 3 mm D x 83 mm H (2 short rails per scale) and 27 in W x .1 in D x 3.3 in H / 686 mm W x 3 mm D x 83 mm H (2 long rails per scale)	
<b>Net Weight</b>	115 lb / 52 kg (including indicator and mounting bracket)	130 lb / 59 kg
<b>Shipping Weight</b>	200 lb / 91 kg (including pallet)	215 lb / 98 kg (including pallet)
<b>Shipping Dimensions</b>	65 in L x 41 in W x 14 in H / 165 cm L x 104 cm W x 36 cm H	
<b>Country of Origin</b>	USA	
<b>UPC Code</b>	809161305005	809161305104

DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.



203 E. Daugherty, Webb City, MO 64870 USA  
Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-2153

**Mobile:** m.detecto.com

**SOLD BY:**